

HIDDEN BLOOD FECAL SEARCH


MATERIAL NEEDED FOR COLLECTION

Specific container provided by the Collection Center (not available for purchase at the pharmacy).



DIRECTIONS AND METHOD OF COLLECTION

1. Collect feces on a clean surface. Avoid contaminating feces with urine.
2. Unscrew and remove the stick with the green cap from the container. (Do not remove the liquid inside the container).
3. Swipe the stick several times over the feces by dipping it in several places. Simply pick up the amount of feces that remains on the grooves of the end of the stick.
4. Insert the stick back inside the container and close it (apply pressure until you hear a "click").
5. Shake the container.
6. Deliver to the Collection Center upon completion of collection. If research is should be performed on three samples, then collect stool on different days and deliver all samples together by storing them in a refrigerated environment.

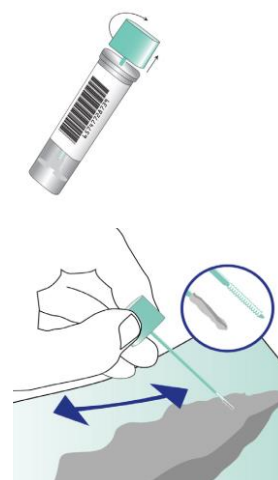

WARNINGS

Do not perform collection during menstruation or if suffering from hemorrhoids.

Warning: container contains sodium azide (< 1%). Keep out of the reach of children. In case of ingestion or contact with eyes and skin, wash the affected area immediately with running water and contact a doctor if necessary.